



---

## Sand Bath 12 lt - S130

### Summary

FALC's sand baths are ideal for the homogenous heating of the contents in the baths. The external structure is made of steel which has been oven coated with acid-resistant epoxy paint whilst the inner tanks are made of extra thick 304 stainless steel.


### Description

FALC's sand baths are ideal for the homogenous heating of the contents in the baths. The external structure is made of steel which has been oven coated with acid-resistant epoxy paint whilst the inner tanks are made of extra thick 304 stainless steel.

The whole structure is insulated with a 40 mm thick layer of ceramic fibre which ensures that the temperature outside the apparatus remains at room temperature.

Heating is provided by iron-clad stainless steel heating elements.

Temperature regulation and control are provided by a bimetallic temperature regulation on a numbered scale. An illuminated green ON/OFF switch is provided as standard.

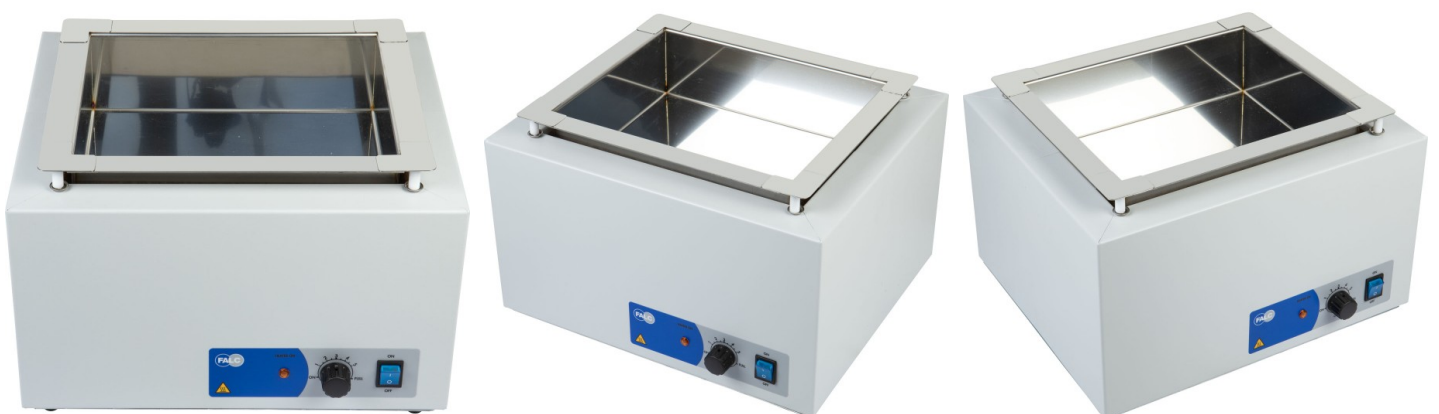

### Features



---

Item	614.1037.13
Temperature	+450°C
Precision	5
Heating power	W 1400
Internal dimension (LxWxH)	mm 500x300x90
External dimension (LxWxH)	mm 610x460x270
Power	V/Hz 230/50-60
Weight	Kg 20.2
Protection class	IP 30
Capacity	lt 12
Sand needed	14 kg